



L-Mesitran®


More than just honey

Dressing guide

Chronic and acute wounds

Superficial wounds

Deep wounds

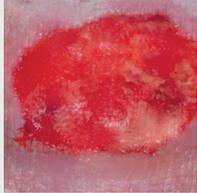

The wound is **RED**

- Protect
- Moist environment
- Bacterial balance
- Stimulate epithelialisation

Clean the wound first (e.g. saline)



- Soft
- Soft
- Hydro / Border
- Ointment
- Soft
- Hydro / Border

- Soft
- Soft
- Soft



The wound is **YELLOW**

- Debride
- Remove bacteria
- Moist environment
- Stimulate healing

Clean the wound first (e.g. saline)

Infected? No

Is the wound wet, moist or dry?




- Soft
- Soft
- Hydro / Border
- Ointment / Soft
- Hydro / Border
- Soft

- Soft
- Soft
- Soft
- Soft

Yes

- Soft

- Soft

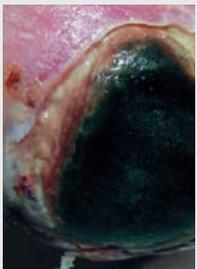

The wound is **BLACK**

- Cleanse
- Remove necrosis
- Remove bacteria

Clean the wound first (e.g. saline)

Infected? No

Is the wound wet, moist or dry?



- Ointment
- Soft
- Ointment
- Soft
- Ointment
- Soft
- Ointment
- Soft

- Ointment
- Soft
- Ointment
- Soft
- Ointment
- Soft
- Ointment
- Soft

Yes

- Ointment
- Soft

- Ointment
- Soft

Post-operative wounds



Closed wounds after suturing

- Ointment
- Soft

Open wounds after surgery

- Ointment
- Soft
- Soft



OINTMENT

- Antibacterial ointment with **48% Medical Grade Honey**.
- In deep or superficial wounds the Ointment will aid debridement and reduce bacterial colonisation. A thin layer is enough for optimal effect.



SOFT

- Antibacterial gel with **40% Medical Grade Honey**.
- Pure honey is known to cause pain when applied to wounds. This is specially designed with sensitive patients in mind due to its lower honey content.



HYDRO

- Hydro-active antibacterial dressing with **30% Medical Grade Honey**.
- The thin honey-hydrogel layer can absorb at least 7 times its own weight of wound fluids. It will follow the contours of the wound and keeps contact with the wound surface. The Hydro provides an ideal healing environment with less pain and reduction of inflammation.



BORDER

- Hydro-active antibacterial dressing with **30% Medical Grade Honey**.
- It's ideal for those difficult locations where regular dressings do not adhere or where mobility is key. It's suited for wounds that need exudate locked away and absorbed, preventing maceration and odour. It will stay on the skin even in the wettest conditions.

Scan me



Visit our website for further information.

L-Mesitran Ointment and Soft are recommended for infected wounds due to the added vitamins and synergistic antimicrobial activity together with honey.